







PULSONIC

SLIM

powered by **BRAUN**

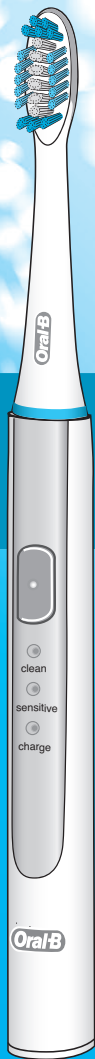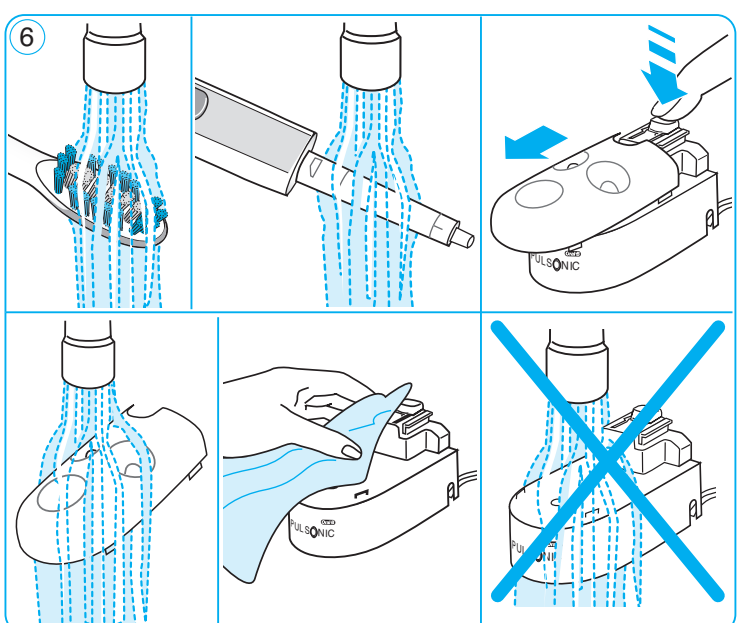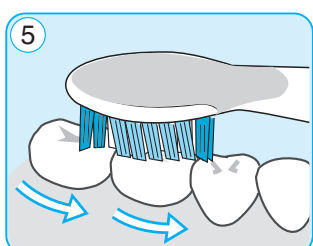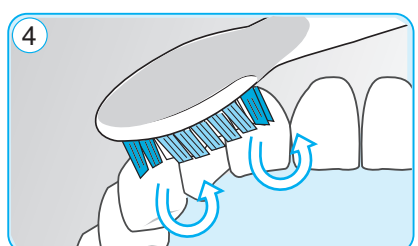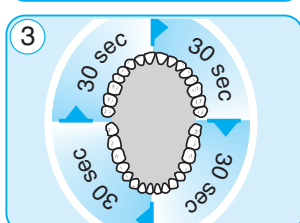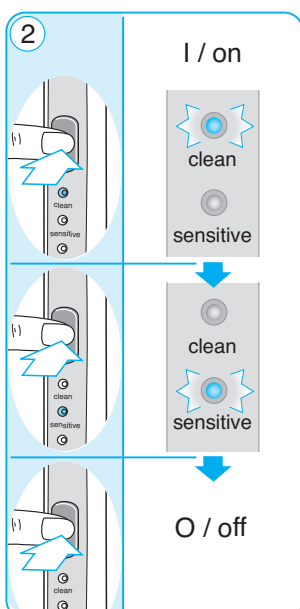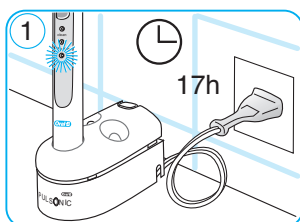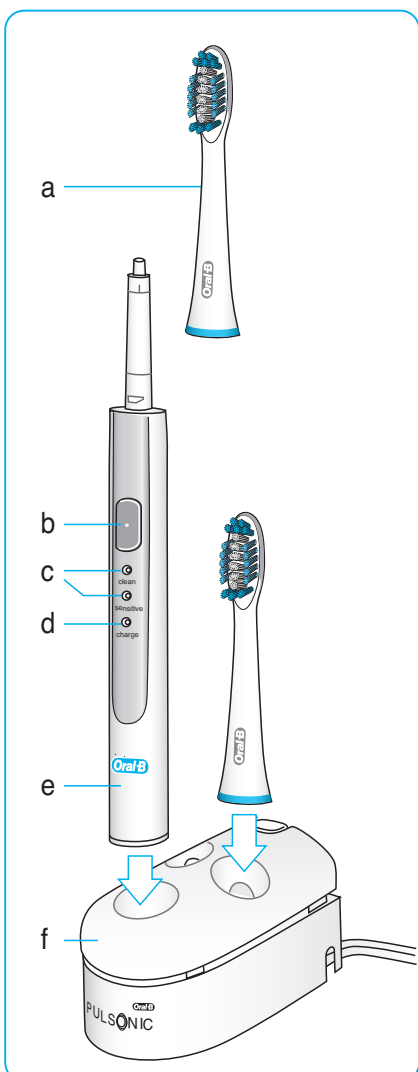Deutsch	3
English	4
Français	5
Español	7
Português	8
Italiano	10
Nederlands	11
Dansk / Norsk	12
Svenska	14
Suomi	16
Ελληνικά	17

Internet:

www.oralb-blendamed.de
www.oralb.com
www.braun.com
www.service.braun.com

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English

Oral-B Pulsonic Slim has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, basic safety precautions should always be followed.

IMPORTANT

- Periodically check the entire product/cord for damage. A damaged or non-functioning unit should no longer be used. If the product/cord is damaged, take it to an Oral-B Service Centre. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water, liquid or store it where it can fall or be pulled into a tub or sink. Do not reach for it if it has fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.

Description

- a Oral-B Pulsonic brush head
- b On/off button
- c Mode indicators
- d Charge indicator
- e Handle
- f Charging unit with handle stand and release button

Connecting and charging

Your toothbrush is electrically safe and designed for use in the bathroom.

- Plug the charging unit (f) into an electrical outlet and place the handle (e) on the charging unit (see picture 1, page 2).
- A full charge takes at least 17 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).

- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power; a battery overcharge is prevented by device. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Using your toothbrush

Brushing modes

Oral-B Pulsonic Slim offers different brushing modes for varying oral care needs:

- «Clean» – Gentle, effective mouth cleaning.
- «Sensitive» – Gentle cleaning for sensitive areas.

When pressing the on/off button (b), your toothbrush automatically starts in the «Clean» mode. To shift from the «Clean» to the «Sensitive» mode, press the on/off button once again. To turn your toothbrush off, press the on/off button (see picture ②, page 2).

Brushing technique

- Slide the brush head on the handle. Wet brush head and apply tooth-paste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on your toothbrush.
- Place the brush head on the outside surface of your upper teeth. The toothbrush bristles should be placed against the teeth at a slight angle towards the gumline.
- Turn your toothbrush on and start brushing in a slightly circular motion (see picture ④, ⑤, page 2). Apply light pressure during brushing. After a few seconds, guide the bristles to the next section. Brush inside, outside and chewing surfaces of your teeth with the same motion throughout the mouth for the recommended 2 minutes.
- To avoid splashing, turn your toothbrush off before removing from the mouth by pressing the on/off button.

You may experience a slight tickling or tingling sensation when you first use your Oral-B Pulsonic Slim. This will subside as your teeth and gums will become accustomed to your Oral-B Pulsonic Slim cleaning experience. If you experience this initial tingling, you may want to brush in the «Sensitive» mode first.

In the first few days of using your toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

Brush head

Oral-B Pulsonic brush head

The Oral-B Pulsonic brush head (a) is recommended for everyday deep cleaning. Its bristles provide a thorough cleaning between your teeth, tooth surfaces and along the gumline. The blue indicator bristles help you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling it is time to replace your brush head.

Timer

A short stuttering sound at 30-second intervals (one break) reminds you to brush all four quadrants of your mouth equally (see picture ③, page 2). A long stuttering sound (two breaks) indicates when you reach the professionally recommended 2-minute mark.

Cleaning

Switch off the handle and remove the brush head. Rinse both brush head and handle thoroughly under running water. Wipe them dry before you put the handle on the charging unit. From time to time you can also clean the charging unit with a damp cloth. For thorough cleaning you may remove the handle stand by pushing the release button on top of the charging unit. The charging unit should never be placed in water (see picture ⑥, page 2).

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste.

For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

30-Days Money Back Guarantee

Try Oral-B starting from the day of purchase.

If you are not satisfied with this product, please return your handle, refills and charging unit in their original box with your original cash register receipt to the address below within 30 days of purchase.

UK

Oral B Products MBG

Dept PX415 , PO Box 13263, Galashiels, TD1 9AJ

ROI

Procter & Gamble

Dept PX415, Oral B Products MBG, PO Box 596, Dublin 14

This does not affect your statutory rights. Terms and conditions apply, for further details, please contact our customer service UK 0800 7311 792 or ROI 1800 409 448 to be referred to the authorised Oral-B Braun Service Centre.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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